

Estd. 1966




DELIVERING TO THE WORLD

Gel & Clot Activator Tube Vacuum (MODEL NO : PPTU 541)



Gel & Clot Activator Tube Vacuum Manufactured by Desco Medical India

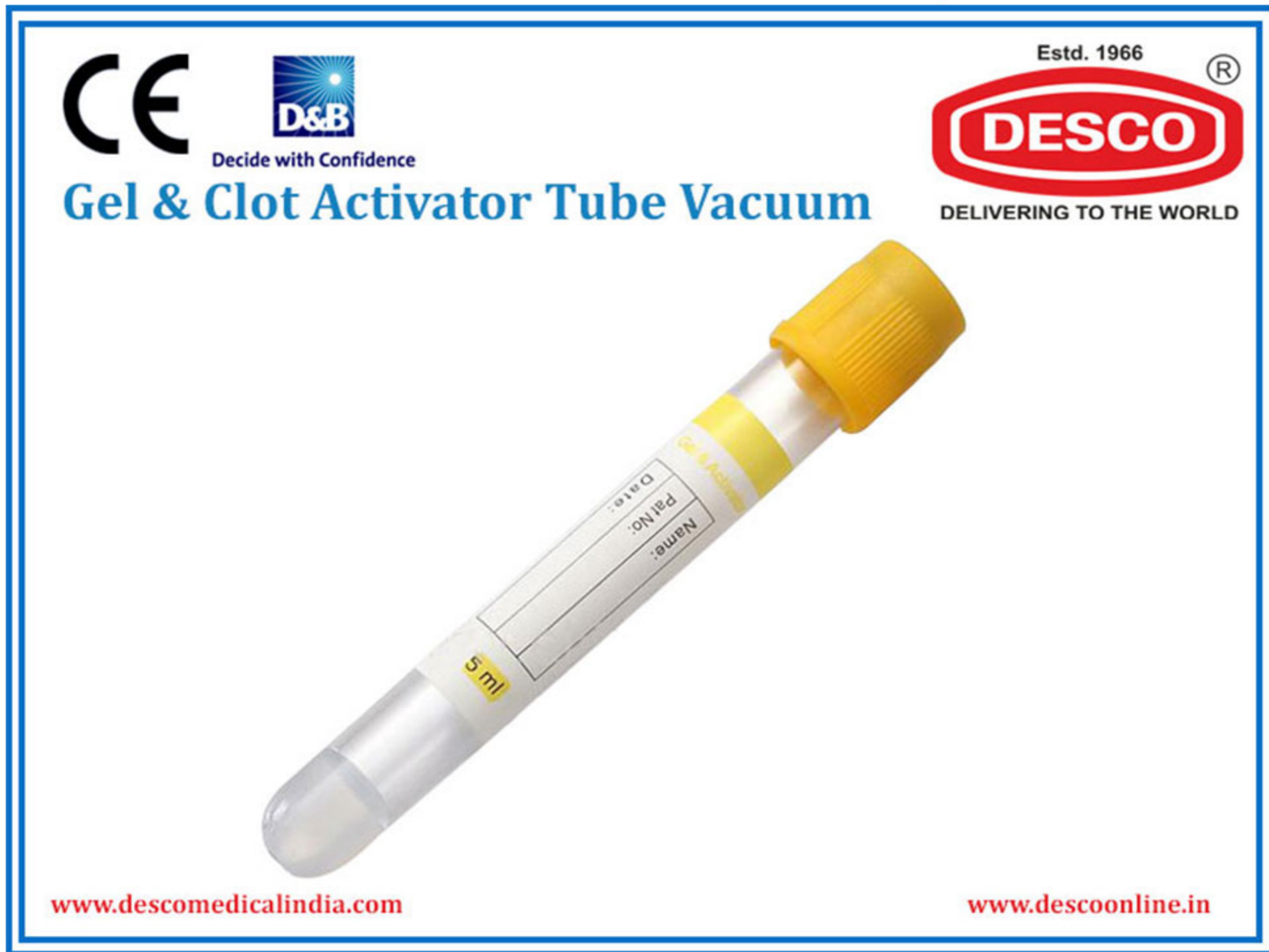

Specifications:-

Gel & Clot Activator Tube is Coated with micronised Silica particles which causes blood to clot rapidly and the barrier gel can effectively separate the serum from fibrin and cells, while preventing substance exchange between blood cell and serum.

Owing to its clot activator, it is better for keeping the original character of the blood and solving the problem of hemolysis.

Cap Colour : Golden Yellow

Max. Vacuum Degree Tolerance : +10 %

Invert the Tube for 6-8 times gently and centrifuge for 10 minutes after 3000-3500rpm after 15-30 minutes coagulation.

It is used in Blood Collection and storage for Biochemistry, Immunology and Serology tests.

We are Manufacturer, Exporters and Suppliers of Gel & Clot Activator Tube Vacuum at wholesale competitive prices from India.

Estd. 1966



DELIVERING TO THE WORLD

Gel & Clot Activator Tube Vacuum (MODEL NO : PPTU 541)



MODEL NO:	OUTER DIA X LENGTH	CAPACITY ML	PACK SIZE
PPTU 541	12.4 x 75 MM	2 ML	100 Pc.
PPTU 542	12.4 x 75 MM	3 ML	100 Pc.
PPTU 543	12.4 x 75 MM	4 ML	100 Pc.
PPTU 544	12.4 x 100 MM	5 ML	100 Pc.
PPTU 545	12.4 x 100 MM	6 ML	100 Pc.
PPTU 546	12.4 x 100 MM	7 ML	100 Pc.
PPTU 547	12.4 x 100 MM	8 ML	100 Pc.
PPTU 548	16 x 100 MM	9 ML	100 Pc.
PPTU 549	16 x 100 MM	10 ML	100 Pc.
PPTU 550	16 x 100 MM	11 ML	100 Pc.
PPTU 551	16 x 100 MM	12 ML	100 Pc.